

IrysView Troubleshooting

Viewing alignments

To view xmaps included in your files move the entire folder that the xmap is in to a Windows machine. Import the xmap you wish to view into IrysView using the import button circled in the screen shot (right). Move all of these files to your windows machine and save them in the same folder but only select the xmap to import and view alignments in IrysView.

Training and installation instructions for IrysView are available here: <http://www.bnxinstall.com/training/docs/>.

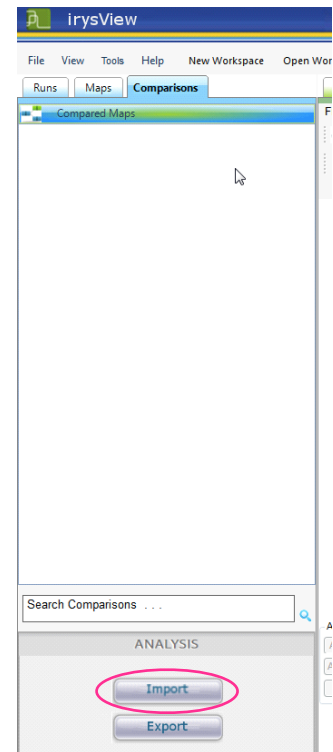


Fig 1: Importing xmaps

Select a reference CMAP to view

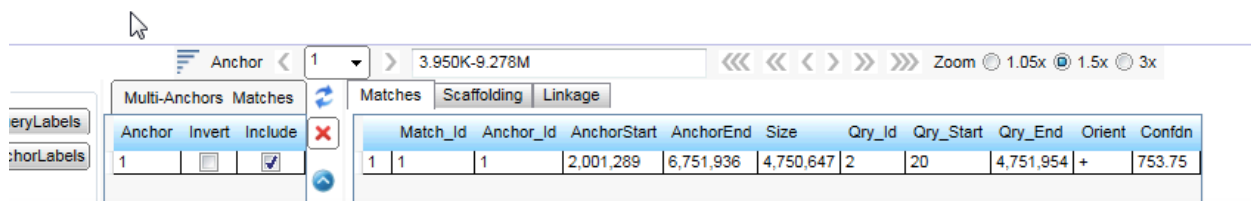


Fig 2: Highlight the XMAP you want to view. Use the “Anchor” dropdown menu to select a single Anchor CMAP (“Anchor” is the same as “Reference”)

IrysView Troubleshooting

If labels don't appear automatically

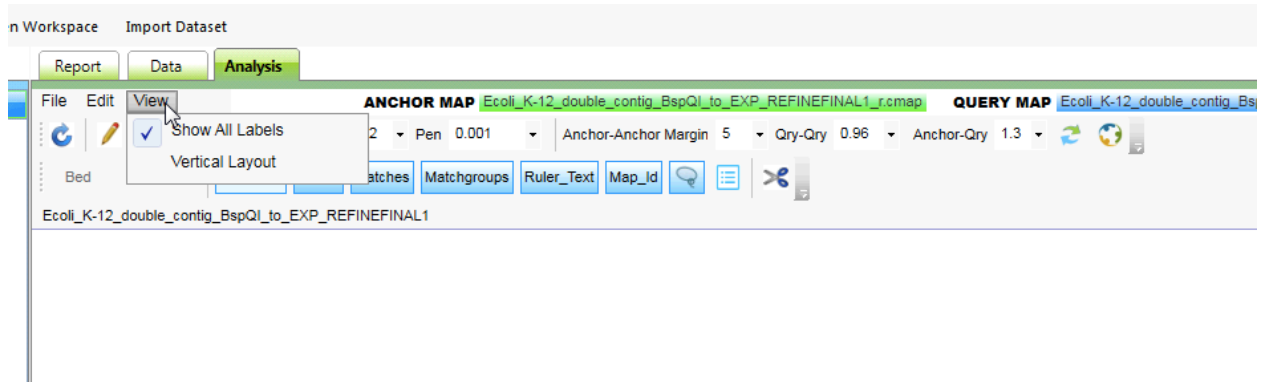


Fig 3: Deselect and reselect “View”->”Show All Labels”. Then press the redraw button with a pencil on it.

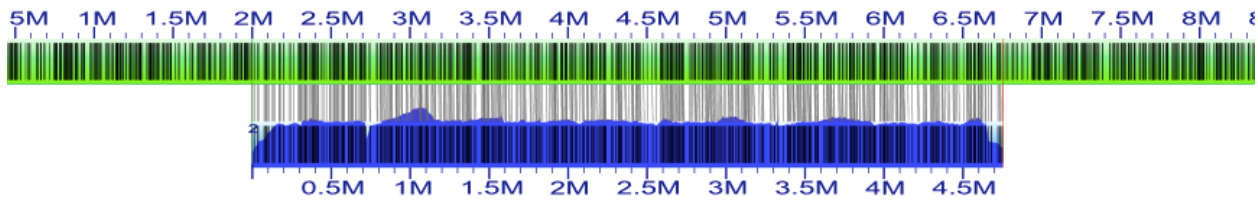


Fig 4: All your labels should now appear.